

**INTERNATIONAL COUNCIL FOR HARMONISATION OF TECHNICAL
REQUIREMENTS FOR PHARMACEUTICALS FOR HUMAN USE**

ICH HARMONISED GUIDELINE

**CONTINUOUS MANUFACTURING OF
DRUG SUBSTANCES AND DRUG PRODUCTS
Q13**

Final version

Adopted on 16 November 2022

This Guideline has been developed by the appropriate ICH Expert Working Group and has been subject to consultation by the regulatory parties, in accordance with the ICH Process. At Step 4 of the Process the final draft is recommended for adoption to the regulatory bodies of ICH regions

Q13
Document History

Code	History	Date
Q13	Endorsement by the Members of the ICH Assembly under <i>Step 2</i> and release for public consultation.	27 July 2021
Q13	Adoption by the Regulatory Members of the ICH Assembly under <i>Step 4</i> .	16 November 2022

Legal notice: *This document is protected by copyright and may, with the exception of the ICH logo, be used, reproduced, incorporated into other works, adapted, modified, translated or distributed under a public license provided that ICH's copyright in the document is acknowledged at all times. In case of any adaption, modification or translation of the document, reasonable steps must be taken to clearly label, demarcate or otherwise identify that changes were made to or based on the original document. Any impression that the adaption, modification or translation of the original document is endorsed or sponsored by the ICH must be avoided.*

The document is provided "as is" without warranty of any kind. In no event shall the ICH or the authors of the original document be liable for any claim, damages or other liability arising from the use of the document.

The above-mentioned permissions do not apply to content supplied by third parties. Therefore, for documents where the copyright vests in a third party, permission for reproduction must be obtained from this copyright holder.

ICH HARMONISED GUIDELINE
CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES AND
DRUG PRODUCTS

Q13

ICH Consensus Guideline

TABLE OF CONTENTS

PART I: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES AND DRUG PRODUCTS.....	1
1. INTRODUCTION	1
1.1. Objective.....	1
1.2. Scope.....	1
2. CM CONCEPTS.....	1
2.1. Different Modes of CM	1
2.2. Batch definition.....	2
3. SCIENTIFIC APPROACHES.....	2
3.1. Control Strategy.....	2
3.2. Changes in Production Output.....	6
3.3. Continuous Process Verification	7
4. REGULATORY CONSIDERATIONS	7
4.1. Description of Manufacturing Process and Process Controls.....	8
4.2. Control Strategy	8
4.3. Batch Description and Batch Size	10
4.4. Process Models	10
4.5. Drug Substance and Drug Product Stability	10
4.6. Conversion of a Batch Process to CM.....	11
4.7. Process Validation	11
4.8. Pharmaceutical Quality System.....	12
4.9. Lifecycle Management	12
4.10. Submission of CM-Specific Information in the CTD.....	12

5. GLOSSARY	15
6. REFERENCES	16
PART II: ANNEXES	18
ANNEX I: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES (CHEMICAL ENTITIES)	18
1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW	18
2. CONTROL STRATEGY AND OTHER TECHNICAL CONSIDERATIONS	19
2.1. Equipment Design and Integration	19
2.2. Process Control and Monitoring	19
2.3. Consideration of Other Controls	21
2.4. Process Validation	21
3. REGULATORY CONSIDERATIONS	22
ANNEX II: CONTINUOUS MANUFACTURING OF DRUG PRODUCTS (CHEMICAL ENTITIES)	23
1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW	23
2. CONTROL STRATEGY AND OTHER TECHNICAL CONSIDERATIONS	23
2.1. Material Characterisation and Control	23
2.2. Equipment Design and Integration	24
2.3. Process Controls and Monitoring	25
2.4. Process Validation	25
3. REGULATORY CONSIDERATIONS	26
ANNEX III: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES (THERAPEUTIC PROTEINS)	27
1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW	27
2. CONTROL STRATEGY	28
2.1. Adventitious Agent Control	28
2.2. Equipment Design and System Integration	28
2.3. Process Monitoring and Real-Time Release Testing	29
3. PROCESS VALIDATION	29
3.1. Approaches to Process Validation	29
3.2. Run Time Considerations	30
3.3. Viral Clearance Validation	30

ANNEX IV: INTEGRATED DRUG SUBSTANCE AND DRUG PRODUCT CONTINUOUS MANUFACTURING.....	31
1. INTRODUCTION	31
2. INTEGRATED SMALL MOLECULE DRUG SUBSTANCE/DRUG PRODUCT PROCESSES.....	31
2.1. Characteristics of Drug Substance and Drug Product Process Steps.....	31
2.2. Example of an Integrated Process.....	31
2.3. Process Design, Monitoring and Control.....	33
2.4. Start-up and Shutdown.....	33
2.5. RTD Characterisation for System Dynamics and Material Traceability	33
3. SPECIFICATION AND BATCH DATA	33
3.1. Drug Substance Specification	33
3.2. Drug Product Specification.....	34
3.3. Example of a Drug Substance and Drug Product Specification	34
3.4. Batch Data.....	35
4. STABILITY REQUIREMENTS.....	35
4.1. Drug Substance Stability	35
4.2. Drug Product Stability	36
5. LOCATION OF DRUG SUBSTANCE AND DRUG PRODUCT INFORMATION IN THE CTD.....	36
ANNEX V: PERSPECTIVES ON MANAGING DISTURBANCES	37
1. INTRODUCTION	37
2. BACKGROUND.....	37
3. MANAGEMENT OF DISTURBANCES	38
3.1. Disturbance Example 1	38
3.2. Disturbance Example 2.....	39
3.3. Disturbance Example 3.....	40

PART I: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES AND DRUG PRODUCTS

1. INTRODUCTION

1.1. Objective

This guideline describes scientific and regulatory considerations for the development, implementation, operation, and lifecycle management of continuous manufacturing (CM). Building on existing ICH Quality guidelines, this guideline provides clarification on CM concepts and describes scientific approaches and regulatory considerations specific to CM of drug substances and drug products.

1.2. Scope

This guideline applies to CM of drug substances and drug products for chemical entities and therapeutic proteins. It is applicable to CM for new products (e.g., new drugs, generic drugs, biosimilars) and the conversion of batch manufacturing to CM for existing products. The principles described in this guideline may also apply to other biological/biotechnological entities.

CM involves the continuous feeding of input materials into, the transformation of in-process materials within, and the concomitant removal of output materials from a manufacturing process. While this description may apply to an individual unit operation (e.g., process chromatography, tableting, perfusion cell culture), this guideline focuses on the integrated aspects of a CM system in which two or more unit operations are directly connected. In this context, any changes made in a unit operation of a CM system may have impact on downstream and upstream unit operations (e.g., back pressure resulting in forward mixing) and output material quality.

Fundamental aspects of CM that are generally not specific to technology, dosage form, or molecule type are described within the main body of this guideline. Annexes are provided to augment the main guideline by providing illustrative examples and considerations specific to certain modalities (e.g., chemical entities, therapeutic proteins), technologies, and production methods (e.g., integration of drug substance and drug product manufacturing). The examples and approaches described in these annexes are illustrative, and alternative approaches can be used. Topics that are broadly applicable to both CM and batch manufacturing are not in the scope of this guideline, and other existing ICH guidelines should be used as appropriate.

2. CM CONCEPTS

2.1. Different Modes of CM

CM can be applied to some or all unit operations in a manufacturing process. Examples of CM modes include:

- A manufacturing approach in which some unit operations operate in a batch mode while others are integrated and operate in a continuous mode

ICH Q13 Guideline

- A manufacturing approach in which all unit operations of a drug substance or drug product manufacturing process are integrated and operate in a continuous mode
- A manufacturing approach in which drug substance and drug product unit operations are integrated across the boundary between drug substance and drug product to form a single CM process (i.e., the drug substance is continuously formed and processed into the drug product through integrated unit operations)

A manufacturing approach may incorporate surge lines or tanks to maintain a constant flow of material inputs and outputs in any mode of CM described above.

2.2. Batch definition

The ICH Q7 definition of a batch is applicable to all modes of CM, for both drug substances and drug products. Based on this definition, the size of a batch produced by CM can be defined in terms of one of the following:

- Quantity of output material
- Quantity of input material
- Run time at a defined mass flow rate

Other approaches to define batch size are possible, if scientifically justified based on the characteristics of the CM process and Good Manufacturing Practice.

A batch size can also be defined as a range. For example, a batch size range can be established by defining a minimum and maximum run time.

3. SCIENTIFIC APPROACHES

3.1. Control Strategy

The development of a control strategy for CM is enabled by a holistic approach, considering aspects specific to CM (discussed below), principles outlined in ICH Q7, Q8, Q10 and Q11, and quality risk management described in ICH Q9.

3.1.1. State of Control

A state of control (ICH Q10) is a condition that provides assurance of continued process performance and product quality. The condition may vary, depending on the mode of CM and the specific process steps. For example, a state of control can be demonstrated for CM processes when a set of parameters (e.g., process parameters, quality attributes) are changing within specified ranges, but the processes are not necessarily in a steady state. Elements of the control strategy monitor the state of control and, when necessary, take appropriate actions to maintain control of the process. It is important to have mechanisms in place to evaluate the consistency of unit operations and the system, and to identify when parameters drift or trend within the specified range. In addition, the root cause of drift or trend, such as variation of inputs, equipment fatigue, or aging of materials, should be identified. For example, in a therapeutic protein process, a change in the elution profile may be due to resin aging.

3.1.2. *Process Dynamics*

Knowledge of process dynamics is important to maintain a state of control in CM. Specifically, understanding the impact of transient events helps to identify risks to material quality and to develop an appropriate control strategy (see Section 3.1.5 for process monitoring and control considerations). Transient events that occur during CM operation may be planned (e.g., process start-up, shutdown and pause) or unplanned (e.g., disturbances) and can be described through time, process parameter and quality attribute values.

Process dynamics should be characterised to understand how output material quality is impacted by transient events. This characterisation could be done by determining properties such as residence time distribution (RTD). RTD characterises the time available for material transport and transformation, and is specific to the process, composition/formulation, material properties, equipment design and configuration, etc. Additionally, understanding RTD and process dynamics enables the tracking of material and supports the development of sampling and diversion strategies, where applicable. Furthermore, such understanding is of importance from a process performance perspective. For example, changes in process dynamics or RTD may impact process characteristics, such as conversion/yield and impurity formation in the manufacture of drug substances.

Process dynamics should be characterised over the planned operating ranges and anticipated input material variability using scientifically justified approaches. Appropriate methodologies (e.g., RTD studies, *in silico* modeling with experimental confirmation) should be used to understand the impact of process dynamics and its variation on material transport and transformation. These methodologies should not interfere with the process dynamics, and the characterisation should be relevant to the commercial process. For example, when conducting RTD studies, the tracer used to replace a constituent of the solid or liquid stream should have highly similar flow properties as those of the constituent replaced, be inert to the other components of the process, and not alter how processed materials interact with equipment surfaces. Step testing by making changes to the quantitative composition of the process stream (e.g., increments of a constituent) is another useful technique to determine the RTD and avoid the addition of an external tracer to the process. Other approaches can be used; the approach taken should be justified.

3.1.3. *Material Characterisation and Control*

Material attributes can impact various aspects of CM operation and performance, such as material feeding, process dynamics, and output material quality. Understanding the impact of material attributes and their variability on process performance and product quality is important for the development of the control strategy. Input materials may require evaluation and control of attributes beyond those typically considered for a material specification used in batch manufacturing. For example:

- In a solid dosage form process, particle size, cohesiveness, adhesiveness, hygroscopicity, static charge, or specific surface area of drug substances and excipients may impact the feeding of powders and material flow through the system.
- In a chemically synthesised drug substance process, viscosity, concentration, or the multiphase nature of the feed may impact flow properties or conversion.

- In a therapeutic protein (e.g., monoclonal antibody) process, lot-to-lot variability of cell culture media or feed components may impact cell culture performance, process performance, or process consistency.

3.1.4. Equipment Design and System Integration

The design of equipment and their integration to form a CM system affect process dynamics, material transport and transformation, output material quality, etc. When developing a CM process and its control strategy, it is important to consider the characteristics of the integrated systems in addition to the individual equipment that can affect process performance. These include the

□ □ □ □ □ □ □ □ □ □ □ □

complete the intended transformation of the material stream within the respective planned operational ranges of the equipment. In addition, transfer steps should also be evaluated, including RTD, for the purposes of integration. Examples of design considerations are given below:

- Design and configuration of equipment (e.g., compatibility and integrity of equipment components for the maximum run time or cycles; geometry of constituent parts to promote the desired transformation; spatial arrangement of equipment to facilitate material flow, maintenance, and avoid build-up or fouling)
- Physical connections and digital control interfaces between equipment (e.g., use of a surge tank between two unit operations to mitigate temporary differences in mass flow rates)
- Locations of material diversion and sampling points (e.g., selection of locations for a diverter valve and sampling probe without disturbing material flow and transformation)

Furthermore, appropriate design or selection of equipment for a CM process may enable process simplification (e.g., through a reduction in the number of unit operations), facilitate process monitoring and material diversion, and improve process capability and performance. For example, in a drug substance process, reactor design can effectively reduce formation and build-up of impurities, resulting in fewer purification steps. Similarly, for therapeutic protein drug substance manufacturing, system design can enable process intensification and reduce cycle times.

3.1.5. Process Monitoring and Control

Process monitoring and control support the maintenance of a state of control during production and allow real-time evaluation of system performance. Common approaches to process monitoring and control including establishment of target setpoints and control limits, design space, and specifications for attributes being measured are applicable to CM.

Process analytical technology (PAT) (ICH Q8) is well-suited for CM. Example applications include in-line UV flow cells to monitor therapeutic protein concentration, in-line near-infrared spectroscopy to assess blend uniformity or water content, and on-line HPLC to monitor conversion of a chemical reaction. The use of PAT enables disturbances to be detected in real time. Therefore, CM is readily amenable to automated process control strategies based on, for example, active process control such as feedforward or feedback process control. Principles of control strategy as described in ICH Q8 and ICH Q11 are applicable to CM processes.

An appropriate sampling strategy is an important aspect of process monitoring and control. The variables monitored, monitoring method and frequency, amount of material sampled (either physical sampling or data sampling using in-line measurement), sampling location, statistical method, and acceptance criteria depend on the intended use of the data and process dynamics. The intended use of data may include detection of rapid changes such as disturbances, assessment of quality of a batch when real-time release testing (RTRT) (ICH Q8) is used, or analysis of process trends or drifts. Further important considerations are the avoidance of measurement interference with the process as well as the impact of physical sampling on the material stream potentially affecting state of control. Assessment of risks associated with data gaps (e.g., PAT recalibration, refill of a feeding system, failure of system components) should inform whether contingency methods are warranted.

3.1.6. Material Traceability and Diversion

CM processes may include periods when non-conforming materials are produced, for example, during system start-up and shutdown or when disturbances are not appropriately managed and mitigated. The ability to detect and divert potential non-conforming material from the product stream during production is an important characteristic of CM and should be considered when developing the control strategy.

Understanding the RTD and process dynamics of individual unit operations and integrated systems over planned operating conditions enables tracking of the distribution of materials over time. This allows input materials to be traced throughout production. Material traceability, understanding how disturbances affect material quality, and the use of appropriate measurements (e.g., PAT) allow for real-time determination of when to start and stop material collection or diversion. The amount of material diverted can be influenced by several factors, such as process dynamics, RTD, control strategy, severity (e.g., magnitude, duration, frequency) of the disturbances, and location of the sampling and diversion points. Additionally, it is important that the diversion strategy accounts for the impact on material flow and process dynamics when material is diverted. Criteria should be established to trigger the start and end of the diversion period and restart of product collection.

3.1.7. Process Models

Process models can be used for development of a CM process or as part of a control strategy for commercial production, including the diversion strategy. Process models may also be used to predict quality attributes in real time, enabling timely process adjustments to maintain a state of control. During development, process models can support the establishment of a design space by explaining how inputs (e.g., process parameters, material attributes) and outputs (e.g., product quality attributes) are related. Process models can enhance process understanding and reduce the number of experimental studies.

For general considerations regarding models (including implications of model impact to validation requirements), refer *Points to Consider: ICH-Endorsed Guide for ICH Q8/Q9/Q10 Implementation*. For CM applications, additional considerations are discussed below.

- A process model is specific to system design and configuration and relevant material properties.

ICH Q13 Guideline

- Model development requires an understanding of the underlying model assumptions (e.g., plug flow versus mixed flow systems) and when these assumptions remain valid. Risk assessments, sound scientific rationales, and relevant data inform the selection of model inputs and model formulation. It is important to determine the relevant inputs that affect the model performance, based on appropriate approaches such as sensitivity analysis.
- Model performance depends on factors such as mathematical constructs and the quality of model inputs (e.g., noise, variability of data). When setting acceptance criteria for model count for uncertainty in the experimental measurement and model prediction should be considered.
- Model validation assesses the fitness of the model for its intended use based on predetermined acceptance criteria using statistically sound approaches. Model validation activities are primarily concerned with demonstrating the appropriateness of the underlying model assumptions and the degree to which sensitivity and uncertainty of the model and the reference methods are understood.
- During commercial manufacture, model maintenance and monitoring of model performance should occur on a routine and ongoing basis considering variability that could impact the model and/or when a process change (e.g., input material, process parameter change) is implemented. Effective and efficient lifecycle management of models is enabled by risk assessment of the impact of a model change (e.g., optimisation of model performance, change of assumptions), considering the scope of model development and model validation criteria. Depending on the extent of a change and its impact on model performance, a model may need to be redeveloped and validated.

3.2. Changes in Production Output

Several approaches to changing production output are discussed below; other approaches are also possible. For already approved products manufactured using a CM process, it is important to justify the selected approach, assess risk and understand its impact on the overall control strategy and process performance, and as needed, update the control strategy to ensure risk control. Some changes may require process modification and process validation.

- **Change in run time with no change to mass flow rates and equipment:** Issues not observed over shorter run times may become visible as run time increases. Additional risks and constraints should be considered and may include, for example, process drift, increased equipment temperatures, material build-up, exceeding the performance limit of components (e.g., validated *in vitro* cell age, resin cycle number), material degradation, membrane or sensor fouling, and microbial contamination. Decreasing production output (below the longest run time previously validated) generally would not imply additional risks, given the same equipment, process and control strategy are used.
- **Increase mass flow rates with no change to overall run time and equipment:** The risks associated with this approach may impact output material quality and are related to changes in process dynamics and system capability to handle increased mass flow rates. Therefore,

this approach may require re-evaluation and modification of the control strategy, including process parameters and controls, material traceability, RTD, sampling, and diversion strategies to reduce risk.

- **Increase output through duplication of equipment (i.e., scale-out):** Considerations for two commonly used scale-out approaches are provided below.
 - *Replication of production lines (like-for-like):* Replicating the integrated CM production line (i.e., same equipment and setup as the original CM system) can be used to increase production output. The replicate production lines follow the same control strategy.
 - *Parallel unit operations on the same production line:* When only some unit operations are replicated on the same line, risks are associated with maintaining control across parallel unit operations. Aspects to consider are maintenance of uniform flow distribution among the parallel operations, synchronisation and re-integration of parallel flow streams, changes to process dynamics, increasing the rate of other unit operations to handle the increased capacity, and material traceability.
- **Scale-up by increasing equipment size/capacity:** Depending on the process and equipment design, increasing production by increasing equipment size may be possible. General principles of equipment scale-up apply, as in the case of batch manufacturing. As elements such as RTD, process dynamics, and system integration may change, various aspects of the control strategy may be impacted. The risks of scale-up and applicability of the original control strategy should be assessed at the new scale and the control strategy should be modified where needed.

3.3. Continuous Process Verification

In CM, frequent process monitoring and control can be achieved through use of process parameters, PAT tools such as in-line/online/at-line monitoring and control, soft sensors, and process models. These tools allow real-time data collection for parameters relevant to process dynamics and material quality, and hence ensure the state of control for every batch. Additionally, since CM can facilitate changes to production output without increasing equipment size, there is an opportunity to generate development knowledge at the same scale intended for commercial manufacturing. These tools, together with the system design and the control strategy, facilitate early execution of process validation activities and the adoption of continuous process verification (ICH Q8) as an alternative approach to process validation.

4. REGULATORY CONSIDERATIONS

The dossier for a CM process should be in accordance with ICH M4Q. Some CM-specific considerations are provided below.

4.1. Description of Manufacturing Process and Process Controls

In line with ICH M4Q, a sequential narrative description of the manufacturing process and process controls should be included in sections 3.2.S.2.2 and 3.2.P.3.3 of the Common Technical Document (CTD) and supported by pharmaceutical development data provided in CTD sections 3.2.S.2.6 or 3.2.P.2. In the case of CM, the information provided in sections 3.2.S.2.2 and 3.2.P.3.3 should be supplemented by the following, when applicable:

- A summary of start-up, shutdown, pause and restart procedures, as applicable, to describe how the integrated CM process functions
- The strategy for material collection and diversion including the criteria that should be met for product collection during routine manufacturing describe

ICH Q13 Guideline

- **Input material attributes:** Impact of input material attributes and their variability (e.g., inter-batch, different suppliers) on continuous processing and product quality should be assessed and proposed material attribute acceptable ranges should be justified when establishing the material specification. For input materials for which pharmacopoeial requirements exist, characterisation and control may extend beyond those requirements.
- **Process monitoring and control:** An appropriate description should be provided in the dossier to show a robust approach to process monitoring and maintaining a state of control. How the control system uses process parameters, in-process material attribute measurements, etc., to make process- and quality-related decisions (e.g., to pause the process or divert material) should be described. Other important aspects should be defined, such as the sampling strategy (e.g., location, sample size, frequency, statistical approach and criteria, and their relevance to the intended use), summary of the models if used (e.g., type of model used, strategies for model evaluation and performance monitoring), and the use of data in making in-process control decisions (e.g., to trigger material diversion).

The data analysis method used should appropriately detect the disturbances or variability in the CM process. For example, when data averaging is used, averaging across appropriate time intervals should be considered rather than data averaging across the entire CM run time. The time intervals should consider the relevancy of PAT measurement frequency to the RTD and process dynamics (i.e., process response time). Therefore, statistical sampling plans and data analysis should be described and justified.

- **System operation:** Procedures should be established and maintained on site for managing system start-up, shutdown, and pauses and for handling disturbances (see Annex V). Relevant approaches for these operations (e.g., handling disturbances) should also be described at an adequate level of detail in the dossier. The disposition of material impacted by transient and pause events should be justified, considering potential risks to output material quality (e.g., the impact of a disturbance as it propagates downstream).
- **Material diversion and collection:** The material diversion and collection strategy should be summarised and justified in the dossier. The strategy described should include the criteria for triggering material diversion, the basis for determining the extent of diverted materials, the conditions for resuming material collection, etc. Factors such as sampling frequency, RTD, and amplitude, duration and propagation of disturbances should be considered in developing the diversion strategy. The extent of diverted material should appropriately incorporate justified safety margins, considering the uncertainty of RTD and other measurements. Procedures for managing material collection, diversion, and disposition (e.g., quarantine, offline testing, investigations) do not need to be included in the dossier but should be maintained within the pharmaceutical quality system (PQS) (ICH Q10).
- **RTRT:** RTRT as described in ICH Q8 can be applied to an output material quality attribute(s). RTRT is not a regulatory requirement for CM implementation. When RTRT is proposed, the associated reference test method should be described. Development of the data collection approach for RTRT implementation should include a risk assessment of

ICH Q13 Guideline

how any lapses in data collection (e.g., recalibrating a near infrared (NIR) probe) may affect decisions relating to product quality. The proposed control strategy should include alternative or additional quality controls to mitigate any risks posed by data lapses. If the results from RTRT fail or are trending towards failure, appropriate investigations should be conducted. Refer to *Points to Consider: ICH-Endorsed Guide for ICH Q8/Q9/Q10 Implementation* for discussion of models used as surrogates for traditional release testing methods.

- **Equipment and system integration:** Aspects of equipment design and system integration that are shown to be critical to output material quality and its control should be described and justified in the context of the overall control strategy.

A summary of the control strategy should be provided in CTD section 3.2.S.2.6 or 3.2.P.2.3 with links or references to the CTD sections that contain the detailed information to enable the understanding and evaluation of the manufacturing process and how it is controlled.

4.3. Batch Description and Batch Size

The approach to define batch size (see examples in Section 2.2) and the proposed commercial batch size or range should be described in the dossier.

If a range is proposed, it should be justified, and the approach for achieving the range should be described (see Section 2.2). Changes in batch size within the approved batch size range can be managed within the PQS. Any post-approval change beyond the approved range should be supported by data (Section 3.2) and appropriately managed (i.e., prior approval or notification).

A suitable quantitative metric with acceptance criteria should be defined within the PQS to establish batch-to-batch consistency and system robustness. For example, when a batch size is defined by the amount of collected material, the amount of diverted materials relative to that of collected materials for each batch should be considered.

The actual intended size of a given batch should be defined before manufacturing begins.

4.4. Process Models

The level of detail provided in the dossier regarding model development, validation, and maintenance over the lifecycle should be commensurate with the model type and impact category. The process model should be specific for the defined system (e.g., equipment, layout, connections). Information to support models used as part of commercial manufacturing should be available during site inspection and maintained at or be accessible to the manufacturing site. Refer to *Points to Consider: ICH-Endorsed Guide for ICH Q8/Q9/Q10 Implementation* for regulatory expectations on process models.

4.5. Drug Substance and Drug Product Stability

Regulatory expectations for the stability data package generally do not differ between CM and batch manufacturing modes (refer to ICH Q1A, ICH Q5C). Batches used to generate primary stability data should be manufactured using a manufacturing process and equipment representative

product and process understanding, system design, and overall control strategy. This justification should be provided in the dossier (see Table 1 below).

When a continuous process verification approach is used to support initial product launch, applicants should justify the appropriateness of the validation activities to provide confidence in the commercial manufacturing process.

4.8. Pharmaceutical Quality System

PQS expectations are the same for batch and CM processes and should follow pertinent ICH guidelines. One important operational aspect of CM is that non-conforming materials can be diverted from the rest of the batch when material traceability, process monitoring, and material diversion strategies are well established. Procedures for material diversion, when required, should be established under the PQS (see Section 4.2). An overarching plan or decision tree that describes how disturbances are managed for various types of material diversion should be maintained under the PQS. Diverted materials resulting from planned events (e.g., system start-up and shutdown) generally do not require investigation when the events meet established process performance criteria. Examples of approaches for managing disturbances are provided in Annex V. As described therein, when unexpected and/or frequent disturbances occur, appropriate investigation, root cause analysis, and corrective action and preventive action (CAPA) should be instituted.

4.9. Lifecycle Management

The principles and approaches described in ICH Q12 are applicable to the lifecycle management of CM.

4.10. Submission of CM-Specific Information in the CTD

The dossier should include information as outlined in ICH M4Q. Additional elements relevant to CM should also be provided in the dossier when applicable. Examples of these elements are listed in Table 1. In the case of integrated drug substance and drug product CM processes, some information and data, such as an integrated flow diagram, may be presented in CTD section 3.2.P with a cross reference in 3.2.S (see Annex IV for additional details).

Table 1: Recommended CM-Related Information in the CTD

CTD section	Information and Data
<p>3.2.S.2.6 3.2.P.2.3</p>	<p>Manufacturing Process Development</p> <ul style="list-style-type: none"> • Summary of the overall process development • Summary and justification of the control strategy, with links or references to the CTD sections that contain comprehensive information on the control strategy elements, for example: <ul style="list-style-type: none"> ○ Strategy for system start-up, shutdown and pause ○ RTD and material traceability information ○ Disturbance management (e.g., material diversion and collection strategy) ○ Sampling strategy, including sampling frequency

ICH Q13 Guideline

	<ul style="list-style-type: none"> ○ Equipment design and system integration aspects that are shown to be critical to output material quality ○ Process controls such as feedforward and feedback process control based on a measured attribute, if used ○ Development and justification of models[#] where used in the control strategy (e.g., process model used for feedforward control, model associated with a PAT analytical method used for in-process testing and/or real-time release) ○ Justification for range of batch sizes, where claimed
<p>3.2.S.2.2 3.2.P.3.2</p>	<p>Batch Definition</p> <ul style="list-style-type: none"> ● Batch size or range, and approach to achieving the intended batch size or range
<p>3.2.S.2.2 3.2.P.3.3</p>	<p>Description of Manufacturing Process and Process Controls</p> <ul style="list-style-type: none"> ● Narrative description of the commercial manufacturing process and flow diagram as described in ICH M4Q, clearly indicating which portion of the manufacturing process is continuous ● Examples of CM-specific aspects of the commercial manufacturing process to be described are: <ul style="list-style-type: none"> ○ Equipment design and system integration aspects when critical to the output material quality ○ CM-related process parameters, controls and criteria (e.g., input rates/mass flow rates, relevant feeder operating ranges), and location points at which process controls or testing is conducted ○ Location of active process controls if used ○ Criteria for product collection, including strategy for diversion
<p>3.2.S 3.2.P</p>	

ICH Q13 Guideline

	critical process parameters and critical quality attributes stay within desired ranges)
3.2.S.4 3.2.P.5	<p>Control of Drug Substance or Drug Product</p> <ul style="list-style-type: none"> • When models are associated with the analytical procedures for release testing of the drug substance or the drug product (e.g., NIR model, dissolution model): <ul style="list-style-type: none"> ○ Summary and justification of the model and the sampling strategy ○ Contingency testing and monitoring plans instituted for when the model is not available (e.g., when gaps in PAT data occur or in case of PAT equipment failure) ○ Analytical model validation information
3.2.R	<p>Regional Information</p> <ul style="list-style-type: none"> • Validation data for high impact process models, if used. • Continuous process verification scheme, if applicable.

The purpose of a model may vary (e.g., testing the quality of in-process material, drug substance intermediate, drug substance or drug product, real-time release testing, process control). Not all categories are covered in this table; information relating to the models should be submitted in the appropriate CTD sections identified in ICH M4Q guideline for these categories. For example:

- Models used for drug substance release testing in 3.2.S.4;
- Models used for in-process testing in 3.2.S.2.4 or 3.2.P.3.4;
- Model used for both in-process testing and real-time release testing in relevant 3.2.S.4 or 3.2.P.5 sections, and incorporated by reference into the applicable Control of Critical Steps and Intermediates section;
- Models used only for process development in 3.2.S.2.6 and 3.2.P.2.3.

5. GLOSSARY

Active Process Controls:

A system consisting of hardware and software architecture, mechanisms, and algorithms that automatically adjust a process to maintain the process output within a desired range. Examples include feedforward and feedback process controls.

Disturbances:

Unplanned changes to process inputs beyond normal operating range or conditions (e.g., process parameter, material property, equipment condition, or environment) that are introduced into a system.

Diversion:

Procedure in which materials are isolated and separated from the product stream in the manufacturing process.

Material Traceability:

The ability to track materials throughout the manufacturing process.

Model Maintenance:

A set of planned activities over the product lifecycle to monitor and performance to continually ensure its suitability for the intended and approved purpose.

Multivariate Statistical Process Control:

The application of multivariate statistical techniques to analyse complex process data with potentially correlated variables. (Ph. Eur.)

Process Dynamics:

The response of a manufacturing process to changing inputs or conditions or transient events.

Residence Time Distribution (RTD):

A measure of the range of residence times experienced by material passing through a specific process environment/vessel/unit operation.

Run Time:

The time interval used to produce a quantity of output material.

Soft Sensors:

A model that is used in lieu of physical measurement to estimate a variable or attribute (e.g., a quality attribute of material) based on measured data (e.g., process data). The model development, including selection of such data variables, is driven by comprehensive product and process understanding.

Steady State:

A stable condition that does not change over time.

System:

A manufacturing architecture that, in the context of CM, consists of individual pieces of equipment, their connections to one another, monitoring and control systems, and spatial layout.

Transient Events:

A temporary condition in which a process goes through a dynamic change. This change may be due to a disturbance or an intentional alteration in the selected operating conditions (e.g., start-up, shutdown, changes from one operating condition to another).

Unit Operation:

A basic step in a process. Unit operations involve a physical, chemical or biological transformation such as: reaction, crystallisation, filtration, blending, granulation, tableting, cell culture, purification or virus inactivation.

6. REFERENCES

ICH M4Q: The Common Technical Document for The Registration of Pharmaceuticals for Human Use: Quality

ICH Q1A: Stability Testing of New Drug Substances and Products

ICH Q5C: Quality of Biotechnological Products: Stability Testing of Biotechnological/ Biological Products

ICH Q5E: Comparability of Biotechnological/Biological Products Subject to Changes in Their Manufacturing Process

ICH Q6A: Specifications: Test Procedures and Acceptance Criteria for New Drug Substances and New Drug Products: Chemical Substances

ICH Q7: Good Manufacturing Practice Guide for Active Pharmaceutical Ingredients

ICH Q8: Pharmaceutical Development

ICH Q9: Quality Risk Management

ICH Q10: Pharmaceutical Quality System

ICH Q13 Guideline

ICH Q11: Development and Manufacture of Drug Substances (Chemical Entities and Biotechnological/Biological Entities)

ICH Q12: Technical and Regulatory Considerations for Pharmaceutical Product Lifecycle Management

Ph. Eur.: European Pharmacopoeia

Points to Consider: ICH-Endorsed Guide for ICH Q8/Q9/Q10 Implementation

PART II: ANNEXES


ANNEX I: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES (CHEMICAL ENTITIES)

1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW

This annex exemplifies one approach to implement CM of drug substances for chemical entities based on the scientific principles described in the main guideline. The examples and approaches in this annex are illustrative, and alternative approaches can be used.

Figure 1 illustrates a drug substance manufacturing process containing both continuous and batch operations. It is not intended to represent a regulatory flow diagram. The continuous process segment consists of unit operations that can be characterised as having two plug-flow reactors (PFRs), two liquid phase extraction operations, continuous crystallisation, and two agitated filter dryers. Manufacture of Intermediate 2 is performed in batch mode, as is final processing including milling and packaging. This annex focuses on the continuous elements of this process.

Figure 1: Example of a drug substance CM system for chemical entities

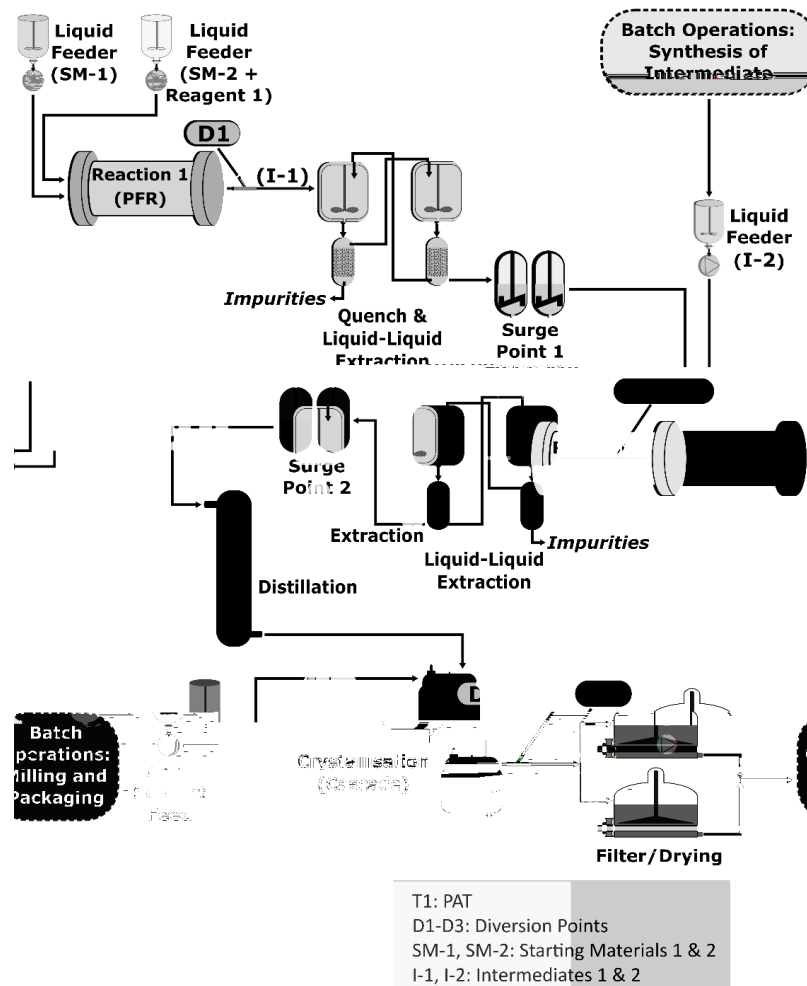

2. CONTROL STRATEGY AND OTHER TECHNICAL CONSIDERATIONS

The CM system and its control strategy were designed to control parameters that impact the manufacture and quality of the drug substance, including impurity profile and physicochemical properties. The overall control strategy was developed in accordance with the main guideline and ICH Q7 Q11.

2.1. Equipment Design and Integration

Within the continuous process segment in Figure 1, the following processes occur:

- Reaction 1: Starting materials 1 and 2 are coupled in a PFR to produce Intermediate 1. Diversion Point D1 is located after the PFR to permit material diversion when PFR conditions are outside predefined acceptance criteria. The reaction is quenched as an integrated operation after the PFR, and unwanted by-products are removed by liquid-liquid extraction. The resultant solution (Intermediate 1) is used as an input for the second reaction without isolation.
- Reaction 2: Intermediate 1 and Intermediate 2 (prepared upstream through separate batch unit operations) are coupled in a second PFR to form the crude drug substance. The online PAT near the reactor exit (T1) monitors conversion of Intermediate 1 to the crude drug substance. Diversion Point D2 located after PAT is used to divert non-conforming material.
- Drug Substance Isolation: The crude drug substance is purified by liquid-liquid extraction and continuous two-stage crystallisation. Distillation prior to crystallisation provides desired concentration of the crude drug substance solution. Diversion Point D3 allows for material diversion at the crystalliser. The crystal slurry is filtered using two agitated filter dryers operating in an alternating fashion to enable continuous operation. The isolated drug substance is then milled using batch operations to achieve the desired particle size distribution prior to packaging.

Two surge points (each containing multiple surge tanks) are used: one before Reaction 2 and another before continuous crystallisation. These are important components of the system design and control strategy, as they improve process robustness and mitigate temporary differences in mass flow rates by decoupling upstream and downstream operations.

The design of the overall system and each unit operation, along with the control strategy, optimise material quality. For example, PFR design elements (i.e., dimension and configuration) allow control of temperature, heat and mass transfer rates, and reaction time. These parameters were shown during development to be important to the drug substance impurity profile.

2.2. Process Control and Monitoring

Holistic controls used across Reactions 1 and 2 ensure consistent operations and quality of the resulting crude drug substance. The stoichiometry of Reaction 1 is controlled via control of concentrations and flow rates of the feeds. Conversion of starting materials to Intermediate 1 with minimal impurity formation is ensured through control of the reaction temperature and flowrate (i.e., reaction time). Reaction 2 is controlled through feedback control of the addition rate of Intermediate 2 based on the PAT measurement of residual Intermediate 1 levels. This ensures

correct stoichiometry for that reaction and minimises the impact of variability of the Intermediate 1 feed solution on drug substance purity. The PAT also measures levels of crude drug substance and impurities, which confirm successful operation of all preceding steps and consistent product quality. Liquid/Liquid extraction ensures appropriate purity control after Reaction 2.

RTD was used to develop a suitable strategy for disturbance detection, corrective actions, and material diversion. RTD characterisation was based on mathematical modeling of all unit operations and surge points across the entire CM process over planned mass flow rates. The RTD was then confirmed through experimental tracer studies for appropriate segments of the commercial equipment. Decisions for triggering material diversion are based on comparing process parameters and PAT measurements to predefined acceptance criteria with timing and duration of diversion informed by the RTD. Importantly, the RTD is also used for material traceability purposes.

Understanding of process dynamics and its impact on quality attributes of material produced throughout the entire process was also used to guide start-up and shutdown strategies. For example, during start-up of Reactions 1 and 2, a small amount of Intermediate 1 or crude drug substance is diverted at Diversion Points 1 or 2, respectively, to allow those materials to reach the target concentrations before processing into subsequent operations. The criteria for diversion were established based on time considering the RTD. This approach was supported by development studies and appropriately confirmed in commercial process equipment. PAT monitoring after Reaction 2 provides additional verification that appropriate criteria have been met during start-up. Collection of material proceeds to the end of the process as subsequently described.

Sampling and process measurement needs were evaluated, considering relevant factors such as residence times (i.e., surge points, process dynamics, and the type and purpose of the measurement). The measurement frequency of the PAT at Reaction 2 enables detection of certain disturbances, informs process adjustments, and ensures timely diversion of material based on predefined criteria. The criteria for material diversion are based on the magnitude and duration of the disturbance, an understanding of RTD and process dynamics of downstream unit operations and surge points, and the impurity purging capability of the crystallisation operation. As a result of this control strategy, all crude drug substance solution that enters continuous crystallisation meets acceptable quality criteria and can be processed through the crystalliser.

Appropriate controls and monitoring requirements for the continuous crystallisation were extensively investigated during development in similar, but smaller scale equipment and appropriately verified using commercial equipment. Process development included spiking studies using impurity-enriched feed solutions and intentional perturbations in process parameters (i.e., feed flow rates, their ratios, and temperatures). An evaluation of the encrusted solids in the crystalliser over extended run times demonstrated the solids were the same form and purity as the free-flowing drug substance slurry. The set of process parameters and ranges identified by these studies were appropriately scaled up. Implementation of these controls along with post-crystallisation material tests (e.g., crystal form, purity) ensure consistent quality of the resulting drug substance throughout continuous crystallisation and subsequent filtration.

The resulting material is milled using a batch operation to provide a drug substance of the appropriate particle size for use in drug product manufacturing. Procedures were developed to allow diversion of material at Diversion Point D3 in the event desired process conditions or material attributes are not met. Diversion of the drug substance from the crystalliser was found to be unnecessary either during start-up or shutdown.

2.3. Consideration of Other Controls

Process robustness and performance over time are important considerations. A risk assessment was performed to ensure that adequate controls are in place to support the proposed run time (which can be up to several months). It identified a number of considerations and corresponding controls/measures. Examples are summarised in Table 2.

Table 2: Examples of other controls for consideration

Consideration	Controls/Measures
Cleaning and fouling potential	<ul style="list-style-type: none"> • Establishment of a risk-based cleaning strategy, including understanding of the impact of build-up on drug substance quality • Additional monitoring to assess fouling and cleanliness (e.g., pressure sensors at the discharge of feed pumps, periodic visual checks for the continuous crystalliser) • Reduction of other risk factors (e.g., filtering feed streams to further reduce fouling risk)
Stability of in-process materials	<ul style="list-style-type: none"> • Hold times at key points in the process (e.g., feed streams; accumulated material at the surge points, reactors, and crystalliser) managed through batch record and process automation • Risk assessment of microbiological growth (i.e., negligible risk based on the nature of the process materials and conditions)
Calibration and potential for changes/drift in instrumentation	<ul style="list-style-type: none"> • Periodic checks at selected points (e.g., process parameter measurements for the PFR, system suitability for the PAT analyser) • Dual sensors at selected locations (e.g., temperature probes for the PFR) so that appropriate corrective actions can be taken
Equipment maintenance	<ul style="list-style-type: none"> • Maintenance requirements for target run time • Use of redundant equipment (e.g., backup pumps) at key locations

Additionally, specifications for input materials were evaluated during process development. There were no differences between batch and continuous processing for this example.

Collectively, the process understanding developed along with implementation of the various controls described provide a robust and reliable control strategy. This ensures consistent quality of the resulting drug substance including the impurity profile, physicochemical properties, and ability of the system to identify and appropriately react to unexpected events.

2.4. Process Validation

The combination of process controls, online PAT measurements, comprehensive monitoring of process parameters and material attributes, and end-product testing results in a data-rich environment for this process. Together with system understanding generated during development,

ICH Q13 Guideline

this enabled the use of a traditional process validation for commercial product launch and continuous process verification to validate process changes over the product lifecycle.

A range of batch sizes was initially established based on material demands and the quantities of material necessary to match input needs of the final batch unit operations. The process was validated using a fixed number of batches. A single planned start-up and shutdown of the commercial CM system was used to manufacture the process validation batches. This approach was supported by the totality of evidence demonstrating the start-up and shutdown capabilities of the system. This included development work on smaller equipment (with appropriate technical justification for being representative), commercial equipment and system qualification data, results of a pre-validation demonstration run, and extensive process monitoring of the CM system that can verify success of each start-up and shutdown in real time.

Subsequently, a


ANNEX II: CONTINUOUS MANUFACTURING OF DRUG PRODUCTS (CHEMICAL ENTITIES)

1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW

This annex exemplifies one approach to implement CM for a tablet drug product based on the scientific principles described in the main guideline. The examples and approaches in this annex are illustrative, and alternative approaches can be used. Specific considerations relating to the implementation of a continuous direct compression process for a chemical entity are presented.

Figure 2 illustrates a continuous direct compression process that consists of continuous feeding, blending, and tablet compression unit operations, with batch mode film coating. It is not intended to represent a regulatory flow diagram.

Figure 2: Example of a tablet drug product CM system

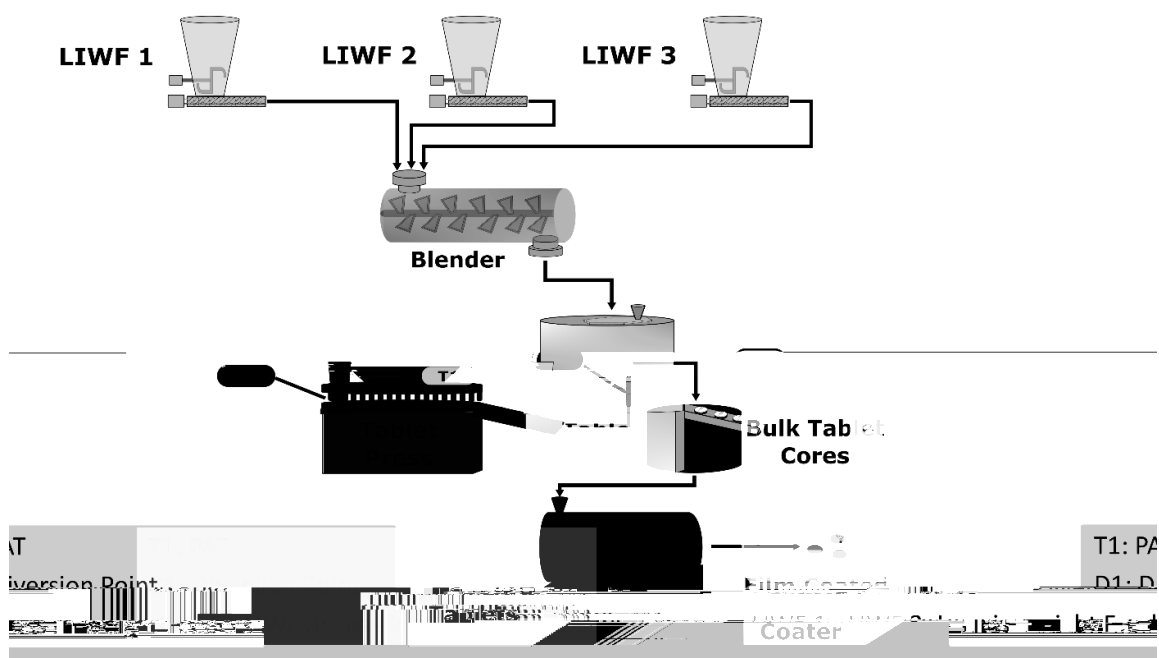

A PAT tool using an NIR method monitors blend uniformity and triggers tablet diversion. Run time at a predefined mass flow rate is used to define the batch size range; in this case, the overall marketing demand requires batch sizes between 360 and 1080 kg of the drug product.

2. CONTROL STRATEGY AND OTHER TECHNICAL CONSIDERATIONS

The CM system and its control strategy were designed to mitigate the impact of disturbances to ensure output material quality. The overall control strategy was developed in accordance with the main guideline and ICH Q8 Q10.

2.1. Material Characterisation and Control

During process design and development, a quality-by-design approach was adopted that identified equipment and process parameters critical to control of the process. Furthermore, the relationships between material quality attributes and their impact on unit operations (particularly the loss-in-weight feeders (LIWFs) and blender) and product critical quality attributes (CQAs) were evaluated. Bulk density of the primary excipient and particle size distribution (PSD) of the drug substance were identified as critical to blend and content uniformity. A defined bulk density range and three-tier (d10, d50, d90) PSD specification were implemented for the excipient and drug substance, respectively.

2.2. Equipment Design and Integration

Unit operations and system components (e.g., NIR probe) were designed or selected to mitigate the impact of disturbances on final product quality. The overall design principle is, where possible, to use gravity to move material. During system integration, the material flow was coordinated across all unit operations to avoid material accumulation or emptying. System mass balance was obtained through understanding of material flow (i.e., RTD) at the intended operating conditions of each unit operation. The impact of equipment design and operation on process dynamics was characterised by the RTD of individual unit operations, as well as the RTD of process segments between individual unit operations and the diversion point. The RTDs were determined by replacing the drug substance in the formulation with a tracer that has highly similar flow properties to those of the drug substance.

The following aspects of equipment design and integration were emphasised:

- **LIWF:** Feeder mass flow rates and their variability were characterised. LIWFs are controlled to deliver the theoretical amount of each input material per the formulation; it was demonstrated that the risks of minor variations to product composition were mitigated by blender mixing capability. Feeder mass flow rates were evaluated using design of experiment (DOE) studies and the proven acceptable ranges of target flow rates were defined. Modelling and statistical approaches were used to help determine the limits for the magnitude and duration of disturbances in mass flow rates, for which material diversion, operator investigation, or process stop are needed. These limits were visualized (e.g., funnel plots) to aid understanding. LIWFs operate in gravimetric mode unless they are refilling (volumetric mode). Refill aspects (e.g., duration and mass of refill) were evaluated to minimise the impact on feeding.
- **Blender:** A horizontal blender was selected for the CM system and the blender design was evaluated (e.g., paddle versus ribbon, number and orientation of paddles in the blender, rotation speed). It was determined that a paddle blender is critical to ensure desired blend uniformity. Rotation speed, number and orientation of the paddles were evaluated for their impact on blend uniformity over the ranges studied, and the corresponding design space for the blending process was defined. RTD characterisation provided information on the degree of forward and back mixing and disturbance propagation, and the RTD was used to define the material traceability and diversion strategy.
- **NIR probe:** The NIR probe was placed in the tablet press feed frame. The chosen NIR equipment met the PAT application requirements (e.g., analysis speed, sampling method,

mass flow rate). Probe location and height are fixed; the impact of material build-up was evaluated and found not significant. The system intended for commercial production was used to generate data for the development, calibration, and validation of the NIR method.

- **Diversion point:** The RTD between the NIR probe and the diversion point was characterised using a tracer. The material diversion strategy links LIWF and NIR limits to the RTD between the LIWF and NIR as well as the RTD between the probe and diversion point, respectively.
- **Coater:** The mass in the coater corresponds to 1 hour of production. Coating was designed to be complete in 45 minutes; whilst coating, the next aliquot of tablet cores is filled into the tablet hopper.

2.3. Process Controls and Monitoring

In this system, the LIWFs may introduce fast dynamic disturbances. These may also occur during changes in operating conditions (e.g., during start-up or process pauses). Therefore, monitoring and control of these events are important elements of the control strategy. The control strategy includes NIR measurements, in-process controls (e.g., individual and total flow rates), process parameters including critical process parameters (e.g., blender rotation speed), and active process controls (e.g., feedback control of tablet weight). The sampling strategy for monitoring and control reflects the observed process dynamics, therefore ensuring adequate detectability of all relevant disturbances. Together, these aspects enable proactive control of the system and ensure continuous operation in a state of control and accurate material diversion to waste based on the predefined criteria. Unique codes are assigned to predefined batch segments to ensure material traceability and identification of conforming and non-conforming materials. Start-up/restart, pause/stop, and shutdown strategies for this example are defined in Table 3.

Table 3: Start-up/restart, pause/stop, and shutdown strategies

Action	Activity
Start-up/Restart	Material tracking and data collection begins; manufactured material is diverted until it meets the predefined acceptance criteria for material collection.
Pause/Stop	A process pause or stop is executed either manually or automatically, according to predefined criteria.
Shutdown	Material collection continues until manufactured material fails the predefined acceptance criteria, and then manufacturing ends.

2.4. Process Validation

In this example, a continuous process verification approach was adopted, considering elements such as prior facility experience in implementing a similar CM process and control system (i.e., platform approach), availability of product-specific data arising from late-stage product development using the commercial equipment, the scale independence of the commercial process (i.e., batch size varies by run time), a comprehensive control strategy with high-frequency data collection, and the use of real-time data from every manufacturing run to further support continuous process verification. The control strategy provides real-time monitoring, trending, and

ICH Q13 Guideline

prediction analysis through the use of NIR measurements, LIWF data, and other data sources arising from monitoring process parameters (e.g., blender torque), thus providing a high degree of assurance of real-time CM system stability and performance and output material quality. The continuous process verification approach, coupled with appropriate regulatory action for reporting manufacturing changes, was used to validate full-time extensions beyond current experience.

3. REGULATORY CONSIDERATIONS

Refer to Section 4 of t-139(d2e 0 2i)-129(4gide)line.9(d2e 0)13(Wn)-139(p)4(ont-9(rsde)r4(a)4(ti)-3(on))-129

ANNEX III: CONTINUOUS MANUFACTURING OF DRUG SUBSTANCES (THERAPEUTIC PROTEINS)

1. INTRODUCTION AND EXAMPLE SYSTEM OVERVIEW


Annex III supplements the main guideline by providing additional regulatory and scientific considerations that are relevant to CM processes for therapeutic protein drug substances and recombinant proteins used as intermediates for subsequent conjugation. It describes aspects that could be applied in fully or partially integrated CM systems. The discussion points and examples presented in the annex are provided for illustrative purposes and are not exhaustive. Alternative approaches can be used.

Figure 3 illustrates a continuous manufacturing process for therapeutic proteins that consists of unit operations characterised by a bioreactor compatible with a perfusion culture system, continuous capture chromatography, virus inactivation, and polishing purification chromatography columns, virus filtration, and buffer exchange and concentration through tangential flow filtration. It is not intended to represent a regulatory flow diagram.

This process integrates a perfusion cell culture with continuous downstream chromatography and other purification steps to continuously capture and purify the target protein. Each individual unit operation is integrated with adjacent unit operations or with a surge line or tank that connects unit operations. Diversion point D1 and PAT (T1) are located after chromatography (Chrom) #1. Using a surge line or tank allows continuous operations to accommodate differences in mass flow rates or process dynamics. Other examples of CM systems may use integrated unit operations for selected steps.

In CM processes, a single thaw of one or multiple vials from the same cell bank may result in either a single harvest or multiple harvests. The number or range of cell bank vials used to produce the specified number of drug substance batches should be defined. The cell bank vials used should be traceable to the output drug substance batches.

Figure 3: Example of a drug substance CM system for therapeutic proteins

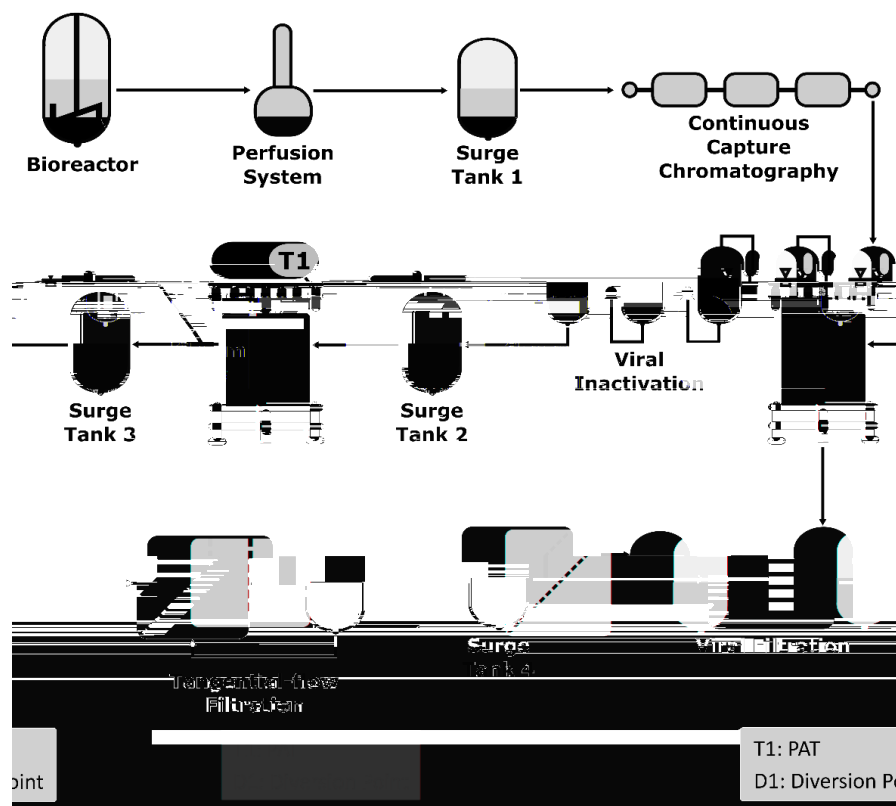

2. CONTROL STRATEGY

2.1. Adventitious Agent Control

In general, all principles used to ensure safety in batch manufacturing are applicable to CM. Safety is demonstrated by a threefold approach based on the principles outlined in ICH Q5A. Control of adventitious agents (e.g., bacteria, viruses, fungi, mycoplasma) should be based on a risk assessment of all potential sources of contamination (e.g., starting and raw materials, manufacturing operations), the ability of the process to remove and inactivate adventitious agents, and the testing capability to ensure the absence of adventitious agents. Based on this assessment, a strategy should be developed to include the type and frequency of adventitious agent testing undertaken to demonstrate that the process remains free of contamination during cell culture and other downstream steps. An aspect unique to CM is extended cell culture duration and continuous processing of harvested cell culture material to obtain drug substances. This means that design measures should be in place to demonstrate the acceptability of all cell culture material used to generate a given drug substance batch. New technologies for real-time decision making, such as rapid testing for adventitious agents, may be used to mitigate the impact of contamination events by earlier detection and appropriate responses during continuous operation.

2.2. Equipment Design and System Integration

The integrity of single-use equipment during use should be ensured to prevent contamination. Single-use connections (e.g., tube welds, connectors) and components may experience long

durations or high change-out frequencies and should be evaluated as potential contamination risks. Filtration steps in CM may be subject to longer filtration periods and potentially increased throughput per unit area or a greater number of filter changes than those under batch manufacturing. Given these factors, a control strategy and a clearly defined scheme should be put in place to allow for filter changes and pre- and post-use integrity testing, as appropriate, without interrupting the process. In the event of a filter failure, a clear strategy for material diversion and refiltration (reprocessing) should be defined.

The CM system should contain appropriate sampling locations based on risk assessment to enable detection of contamination, while avoiding unnecessary contamination risk introduced through the sampling procedure. The sampling locations and frequency may be adjusted based on improved product and process understanding.

Integrated systems may use surge tanks for flow rate adjustments or other purposes between processing steps. When surge tanks are used, the relevant RTD, uniformity and microbial risks to the product in these surge tanks should be evaluated and defined in advance.

2.3. Process Monitoring and Real-Time Release Testing

CM lends itself to various monitoring schemes with different levels of automation. Examples include in-line sensors placed directly in a process vessel or flowing material stream and online analysers that conduct automatic sampling. Regardless of the approach used, appropriate monitoring at suitable stages of the CM process enables timely data analysis to ensure operations are in a state of control. In certain cases, relevant process parameters may be adjusted to ensure the quality of in-process or output materials. Enhancing in-line/online PAT capabilities and development of automation systems for process monitoring enable a continuous monitoring scheme in support of a release testing strategy that may include RTRT for some quality attributes. For example, drug substance in-line release tests for pH, osmolality, protein concentration and online release tests for purity, charge heterogeneity, aggregation, and low-molecular weight impurities can be performed at specific points in the drug substance manufacturing process shown to be critical for control of the product quality attributes.

An in-process test that is demonstrated to be representative of a release test performed on the output material may be used as a surrogate for a traditional release testing method. Conventional offline testing for product release is necessary for quality attributes for which analytical technologies are not available for online or in-line measurements (e.g., potency). Likewise, conventional tests for monitoring and control (e.g., microbiological analytical methods and other tests that require long processing times) might also be needed.

3. PROCESS VALIDATION

3.1. Approaches to Process Validation

Process validation approaches used for batch processes are also applicable to CM processes with additional considerations for equipment performance over extended run times, integration of automation systems, and material flow throughout a CM system. Therefore, the scope of validation continues to demonstrate the ability of the CM system to consistently manufacture a product with the desired quality.

For therapeutic protein CM, any approach chosen to demonstrate the consistency of process performance and product quality should also consider sources of variability that may have a potential impact on product quality. This may include variability between batches purified from harvest materials collected up to the limit of *in vitro* cell age from a single cell bank thaw, as well as the potential variability between different batches purified from harvests of multiple cell bank thaws. Variability may be evaluated either as part of process validation or through alternative studies, if justified.

Alternative process validation approaches (e.g., continuous process verification) may be considered when justified. Elements such as risk assessment, applicability of development data, control strategy, and prior knowledge may be considered in determining the suitability of an alternative process validation approach.

3.2. Run Time Considerations

Bioreactors for CM may operate for significantly longer periods of time than bioreactors for batch manufacturing. The approach to establish a limit of *in vitro* cell age for production cells does not differ, regardless of the mode of bioreactor operation. Previously established limits of *in vitro* cell age for a bioreactor operating in a batch mode run may not be applicable to a bioreactor operating in a continuous mode under different culture conditions. The limit of *in vitro* cell age used for production should be based on data derived from production cells expanded under pilot-plant scale or commercial-scale conditions to the proposed *in vitro* cell age or beyond as outlined in ICH Q5A, Q5B and Q5D.

Run time considerations should include factors such as the control of all adventitious agents (e.g., viruses, bacteria, fungi, mycoplasma) and resin and membrane lifetimes. Viral testing should be conducted as outlined by ICH Q5A, and an appropriate microbial control strategy should be established.

3.3. Viral Clearance Validation

The general recommendations outlined in ICH Q5A for viral safety and clearance remain applicable for CM. Where recommendations may not be applicable to a CM system, scientifically justified alternatives may be proposed.

ANNEX IV: INTEGRATED DRUG SUBSTANCE AND DRUG PRODUCT CONTINUOUS MANUFACTURING

1. INTRODUCTION

Annex IV supplements the main guideline by providing additional regulatory and scientific considerations that are relevant for the development and implementation of integrated drug substance and drug product CM processes (referred to as integrated process(es) hereafter).

This annex also provides an example of an integrated process for a small molecule tablet dosage form. The example and approaches described in this annex are not exhaustive. Alternative approaches can be used.

2. INTEGRATED SMALL MOLECULE DRUG SUBSTANCE/DRUG PRODUCT PROCESSES

2.1. Characteristics of Drug Substance and Drug Product Process Steps

Considering the differences between the drug substance and drug product process steps enables appropriate design of an integrated process. For example, process steps for drug substance and drug product manufacturing can have different RTDs, and a prevalence for liquid or solid input material addition can lead to a different frequency of in-process measurements. These differences are expected to influence the selection of equipment, equipment connections, surge lines or tanks, and the locations of in-process measurements and material diversion.

2.2. Example of an Integrated Process


Figure 4, which is not intended to represent a regulatory flow diagram, illustrates a fully continuous integrated drug substance and drug product process. It shows the following elements:

- Material addition points for liquids and solids
- Each process step used for drug substance and drug product manufacturing
- Process design for the interface between the drug substance and drug product
- Sampling locations for all in-line/at-line/offline measurements, including PAT (shown by T1 T5)
- All diversion points (shown by D1 D4)

In this example, chemical reaction using flow reactors, continuous crystallisation, and crossflow filtration are used to obtain the drug substance as a highly concentrated crystal slurry. A wet granulation process consisting of blending, granulation, drying, milling, compression and coating unit operations is used to obtain a tablet drug product. The selection of a wet granulation process for the manufacture of the drug product permits the drug substance and drug product processes to be integrated through the continuous filtration line. The concentrated crystal slurry functions as both the drug substance source and the granulation fluid. No surge lines or tanks are used.

Other process schemes including, for example, different purification methods, surge tanks, mix of batch and continuous unit operations could also be used in the design of an integrated process. Details should be provided on how drug substance purity is ensured.

Figure 4: Example of an integrated drug substance and drug product CM system

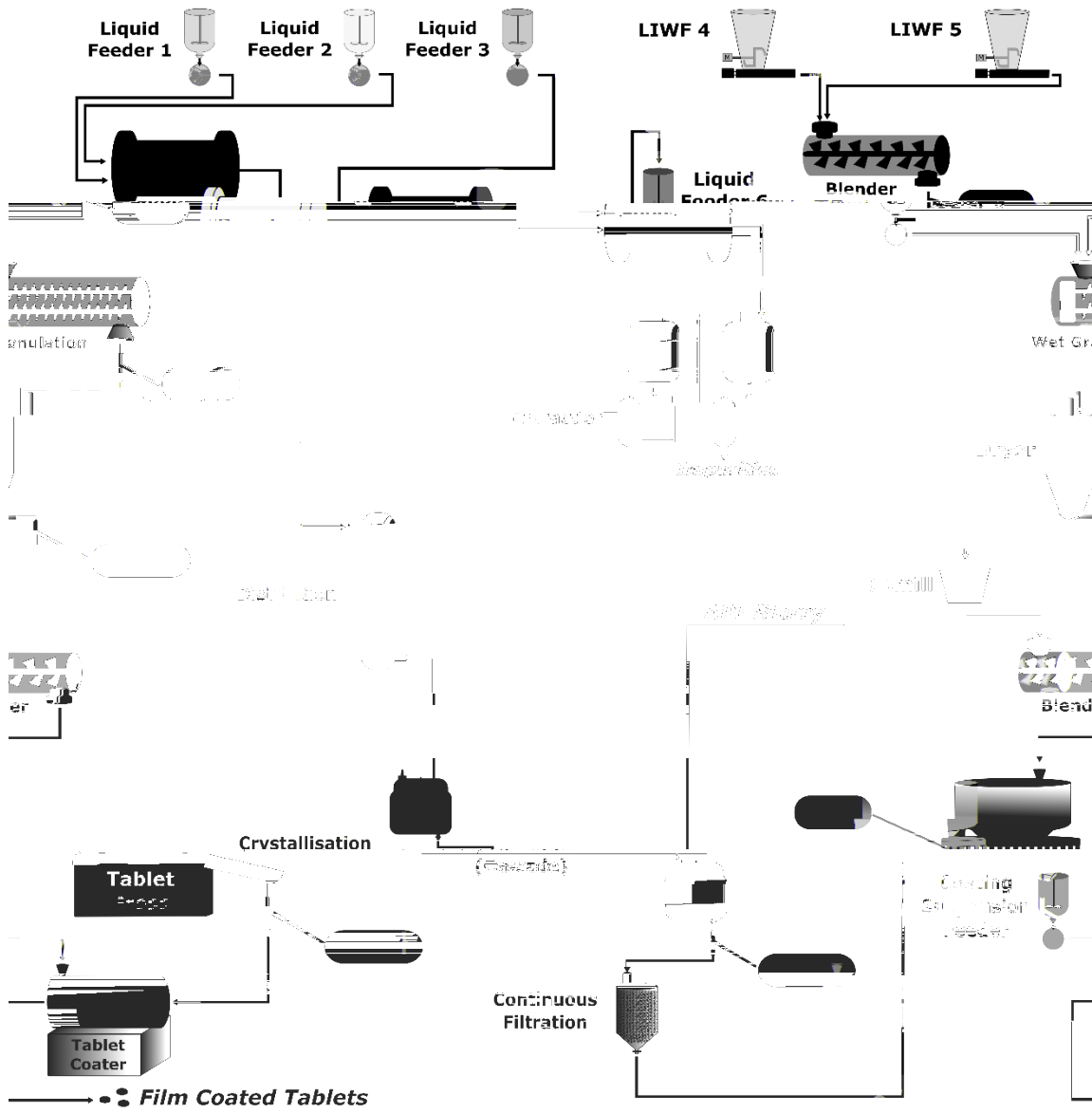

Locations
Liquid Feeders

T1 - T5: PAT and at-line test loc
D1 - D4: Diversion Points
LIWF 4, LIWF 5, have in weight

facility changes), investigation of adverse events and product recalls, development of pharmacopeial monographs, and the establishment of a reference standard.

Although a drug substance specification should be instituted, drug substance testing may not be required on a routine basis when the integrated process is appropriately controlled. A set of process performance criteria comprehensive monitoring of the quality of the drug substance during the lifecycle of the product, conformance to the drug substance specification should be verified on a periodic and event-driven basis by testing the purified drug substance at an appropriate location using a relevant sampling plan. The frequency of the periodic verification should be defined and justified. Drug substance periodic verification can be based on pre-selected batches and/or at predetermined time intervals (e.g., every 5th batch or if manufactured infrequently, annually). Event-based verifications could be triggered by a change in supplier, starting material, synthesis conditions, or other factors considering risk. Refer to ICH Q6A for additional details on periodic testing.

Appropriate sampling locations should be incorporated into the process design to enable testing of the drug substance (e.g., location T1 in Figure 4). Any modifications made to the sample to enable the test (e.g., drying of the crystal slurry for testing crystalline form) should be incorporated into the test methodology. Sampling locations should be identified in the drug substance specification.

Although the drug substance is not isolated, the origin and fate of potential impurities (e.g., related substances, residual solvents, catalysts), robustness of impurity clearance, and impurity carryover from the drug substance into the drug product should be provided in the dossier. The control of impurities formation and clearance should be integrated into the overall control strategy.

3.2. Drug Product Specification

In integrated processes, attributes typically associated with the drug substance quality are generally included in the drug product specification unless justified per ICH Q6A. Therefore, the drug product specification in an integrated process should include drug substance-related substances, residual solvents (used in drug substance synthesis), elemental impurities, etc., when appropriate. The specified impurities in the drug product specification may differ from the specified impurities in the drug substance specification.

Sampling location should be appropriately identified in the drug product specification table, as some testing (e.g., testing for drug substance periodic verification as described above) may need to be performed following the drug substance purification step (before drug product formation).

3.3. Example of a Drug Substance and Drug Product Specification

An example of a drug substance and drug product testing approach for an integrated process based on Figure 4 is shown in Table 4. The test attributes listed are considered relevant for this example. The specific details of each integrated process should be considered in the selection of the appropriate test attributes and testing plan.

Table 4: Example of a testing approach for an integrated CM

Test Attribute ¹	Drug Substance Specification for Periodic Testing		Drug Product Specification for Routine Testing of Every Batch	
	Test	Sampling Location	Test	Sampling Location
Description	N/A	N/A	✓	Coated tablets
Identity	✓	Use drug product test result	✓	PAT at tablet feed frame (T4)
Crystalline Form ²	✓	Sampling Location T1	N/A	Not tested when justified
Chirality ³	✓	Sampling Location T1	N/A	Not tested when justified
Particle Size	✓	Sampling Location T1	N/A	Not tested
Purity	✓	Sampling Location T1	N/A	Not tested
Assay	N/A	N/A	✓	Core tablets, Sampling Location combination of T4 (blend uniformity) and T5 (tablet weight)
Impurities ⁴	<i>Impurity specification for drug substances and drug products may differ</i>			
<i>Related Substance</i>	✓	Sampling Location T1 (at-line high performance liquid chromatography (HPLC)) ⁴	✓	Sampling Location T1 (at-line HPLC) ⁴ or Coated Tablets (offline HPLC testing), as appropriate
<i>Residual Solvents</i>	✓		✓	
<i>Elemental Impurities</i>	✓		✓	
<i>Mutagenic Impurities</i>	✓		✓	
Dissolution	N/A	N/A	✓	Coated Tablets
Uniformity of dosage units	N/A	N/A	✓	Uncoated Tablets
Water content	N/A	N/A	✓	Coated Tablets
Microbial limits	N/A	N/A	✓	Coated Tablets

for other aspects, such as to support the storage of in-house reference standards and to gain an understanding of product stability profiles. Institution of a hold time enables temporary storage of drug substance during an interruption in production.

4.2. Drug Product Stability

The ICH stability guidelines and Section 4.5 of the main guideline are applicable to drug product made by an integrated process.

5. LOCATION OF DRUG SUBSTANCE AND DRUG PRODUCT INFORMATION IN THE CTD

Drug substance and drug product information could be provided in the respective CTD sections 3.2.S and 3.2.P of the dossier. A description of the process step that integrates the drug substance and drug product could be based on its relevancy to the respective section. For example, in the process example provided in this annex, the continuous filtration process could be described in CTD section 3.2.S as it is related to concentration of the drug substance. The integrated flow diagram can be provided in CTD section 3.2.P and referenced in section 3.2.S.

ANNEX V: PERSPECTIVES ON MANAGING DISTURBANCES

1. INTRODUCTION


This annex describes examples of approaches for managing transient disturbances (hereafter referred to as disturbances in this annex) that may occur during CM. The discussion points presented here are not exhaustive. Alternative approaches can be used.

2. BACKGROUND

Disturbances may result in a variation in material quality. The impact of some variations on material quality in an earlier process step may be resolved by downstream process steps. The extent of variations and the ability to resolve them in subsequent steps are impacted by the amplitude, duration, and frequency of the disturbance. Identification of tolerable ranges for these parameters and establishing appropriate acceptance criteria will enable the development of an effective strategy for managing disturbances.

Manufacturers may use various methodologies (e.g., DOE, RTD studies, a combination of both) to understand the impact of disturbances. Funnel plot predictions based on an RTD model can be a useful tool to understand the qualitative and quantitative impact of the amplitude and duration of a disturbance on material quality. Figure 5 shows a funnel plot for drug substance feeding in a drug product CM process (similar to the example in Annex II). Funnel plots are specific to the formulation, process conditions, and system configuration used in RTD model development. Information from the funnel plot helps to inform the selection of appropriate acceptance criteria for disturbances. For example, the dotted lines in the following funnel plot show that a disturbance of +/- 20% lasting less than 90 seconds would not cause the drug concentration in the blend to exceed the 90 110% label claim (LC).

Figure 5: Example of a funnel plot for the feeding of a drug substance

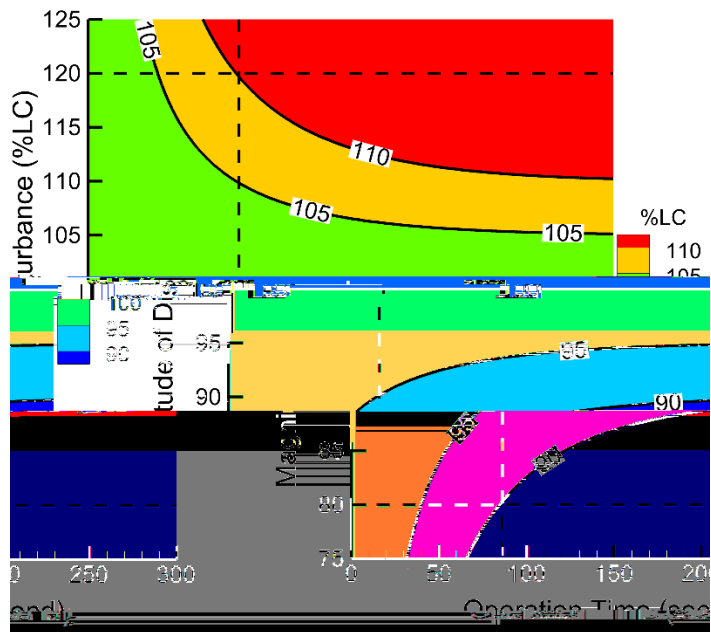

3. MANAGEMENT OF DISTURBANCES


Manufacturers may develop various approaches to manage disturbances considering the specific details of the CM system and the risk to material quality from disturbances. Three examples considering different risks of a disturbance are provided below:

- Example 1: The amplitude and duration of the disturbance meet predefined acceptance criteria for the disturbance, and the occurrence of such disturbances is infrequent.
- Example 2: The amplitude or duration of the disturbance exceeds the predefined acceptance criteria for the disturbance, and the occurrence of such disturbances is infrequent.
- Example 3: The amplitude and duration of each disturbance meet predefined acceptance criteria for the disturbance, but multiple, frequent disturbances are observed.

These examples focus on the impact of disturbance from a loss-in-weight feeder (LIWF) on the drug concentration in the blend for a CM process similar to that described in Annex II, given that all other parameters being monitored meet the predefined acceptance criteria. These examples use the information in the funnel plot (Figure 5) and, for the purpose of discussion, assume that the acceptance criteria for the magnitude and duration of an LIWF disturbance is +/- 20% lasting for 80 seconds. These examples help illustrate the important considerations in management of disturbances under selected scenarios, which may also be applicable to drug substances and other CM processes.

3.1. Disturbance Example 1

Figure 6: Example of an infrequent disturbance that is within the acceptance criteria for disturbances

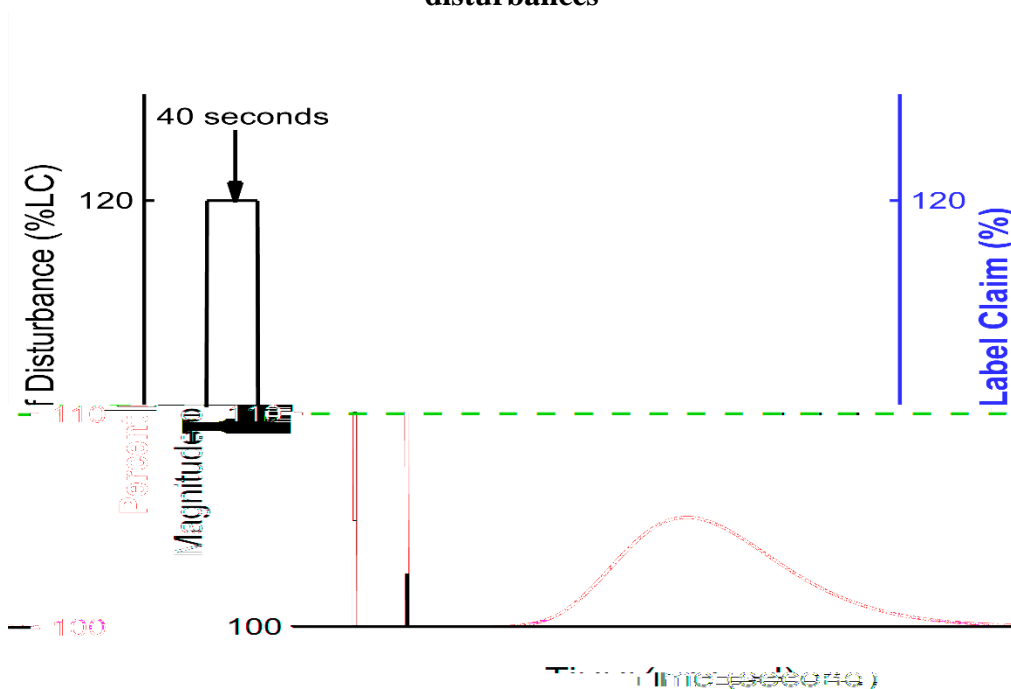


Description: Figure 6 illustrates a drug substance LIWF with an infrequent transient +20% flow spike lasting 40 seconds, which is within the predefined acceptance criteria for disturbances. This disturbance causes an increase in the amount of the drug substance fed into the blender, before returning to normal operating condition. The funnel plot (Figure 5) shows that following this disturbance, the drug substance concentration in the blend remains within the 90-110% acceptance criteria, due to back mixing. An additional quality check, such as measurement of the drug substance concentration at a suitable location (e.g., NIR measurements at the tablet press feed frame) could be considered, to confirm the blend is within 90-110%.

Impact: Although this disturbance represents an excursion from normal operation, the quality of the output material is not affected as the magnitude/amplitude of the disturbance and product quality meet their predefined acceptance criteria.

Action: No material is diverted. Collection of the output material continues, and the process continues to operate. No investigation is needed, because such a disturbance has been evaluated during development and impact on material quality is understood.

3.2. Disturbance Example 2

Figure 7: Example of an unexpected disturbance that is outside the acceptance criteria for disturbances

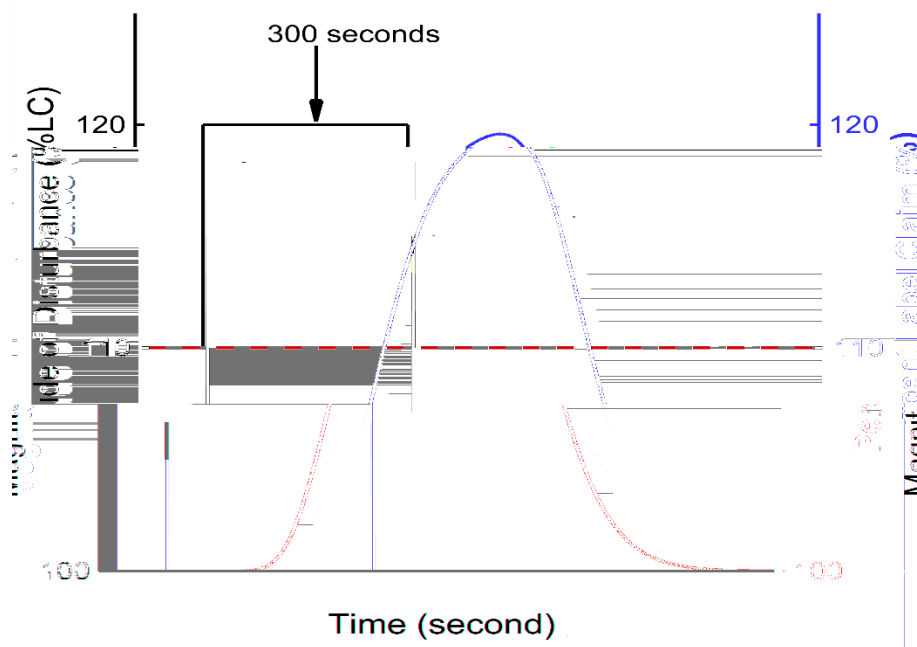

Description: Figure 7 illustrates a drug substance LIWF with an unexpected transient +20% flow spike lasting 300 seconds. The disturbance is outside the predefined acceptance criteria for disturbances. This disturbance causes an increase in the amount of the drug substance fed into the blender before returning to normal operating condition. The funnel plot (Figure 5) shows that following this disturbance, the drug substance concentration in the blend exceeds the 90-110%

acceptance criterion. An additional quality check, such as measurement of the drug substance concentration at a suitable location (e.g., NIR measurements at the tablet press feed frame), confirms the blend exceeds 110%.


Impact: The quality of the output material is adversely impacted as the disturbance duration exceeds the predefined acceptance criteria.

Action: The process continues to operate while the non-conforming material is diverted according to a pre-established procedure, and the time to start and end diversion is controlled by the automation system. The system returns to normal material collection mode when the non-conforming material is completely diverted. If needed, a concurrent investigation is initiated to determine root cause.

Diverted Amount: The amount of material diverted depends on the control strategy used (including specific triggers for material diversion) and on the process dynamics from the point of disturbance detection and the point at which material diversion ends. Inclusion of confidence intervals in the RTD provides a safety margin to ensure all non-conforming material is diverted from the batch. Additional factors, such as the sampling strategy and the ability to trace and remove materials, are considered in establishing the criteria for material diversion.

3.3. Disturbance Example 3

Figure 8: Example of disturbances that are within the acceptance criteria for disturbances, but occur frequently

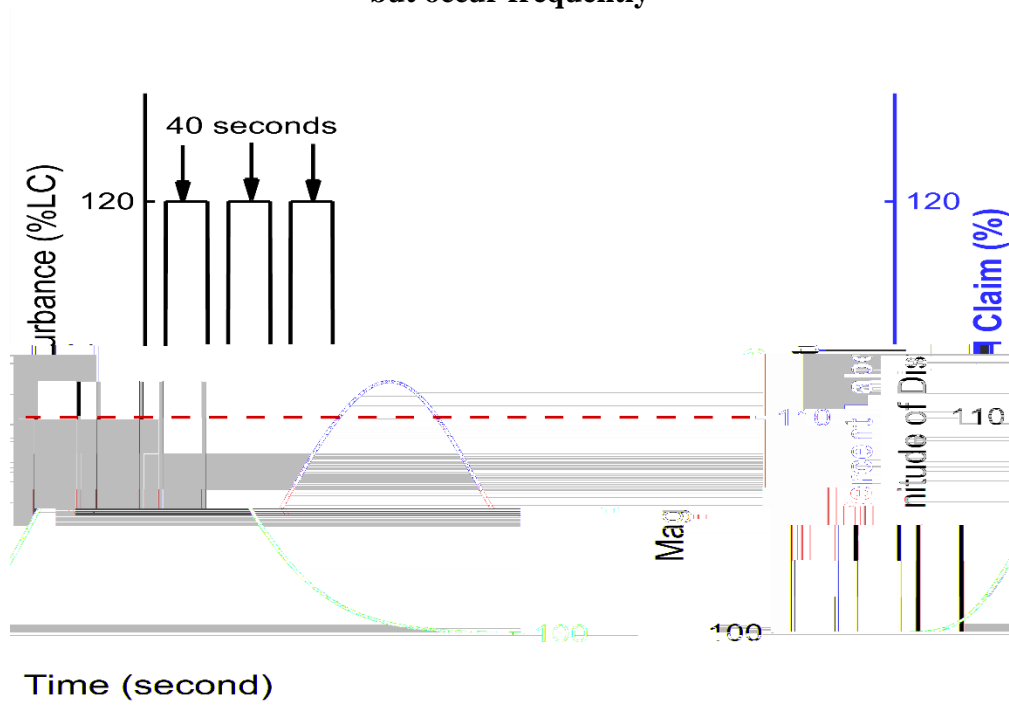

Description: Figure 8 illustrates a drug substance LIWF with multiple frequent transient +20% flow spikes, each lasting 40 seconds, resulting in variability in the amount of material fed into the blender.

ICH Q13 Guideline

Impact: Although each disturbance meets the predefined acceptance criteria for disturbances, it occurs with a high frequency over a short time period. In this example, the system cannot dampen these multiple disturbances sufficiently, thus resulting in non-conforming materials.

Action: The impact of these disturbances on system performance and output material quality is monitored closely (e.g., NIR method, other elements of the control strategy). Process operation and product collection continue until one or more elements of the control strategy do not meet the predefined acceptance criteria. When a criterion is no longer met, the material is diverted according to a pre-established procedure. If high-frequency disturbances persist, process operation may be paused. An investigation is conducted to understand the root cause for these frequent disturbances. Such investigations enable preventative actions to be taken to avoid equipment failure and adverse impact on critical quality attributes, ensure process performance (e.g., robustness), etc. Assessment of process capability or other evaluations may also be warranted. Setting acceptance criteria for the frequency of disturbances could also be considered to aid the management of disturbances.

Diverted Amount: The amount diverted is the same as described in Section 3.2 of this annex. The disposition of the diverted material and the entire batch is assessed upon completion of the investigation.